




# **To study electro-deposition of Cu on cathode**

# Electro-deposition or Electroplating

- The chemical conversion of ions in solution into a solid deposit of metal atoms with the work of a electrical power supply

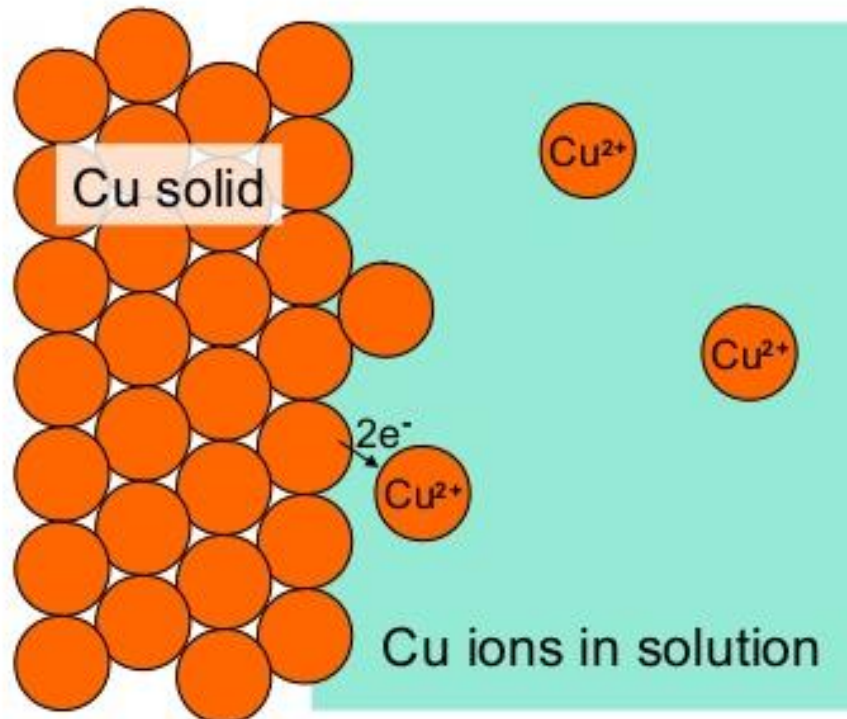


# Electro-deposition cell – Electrolytic cell

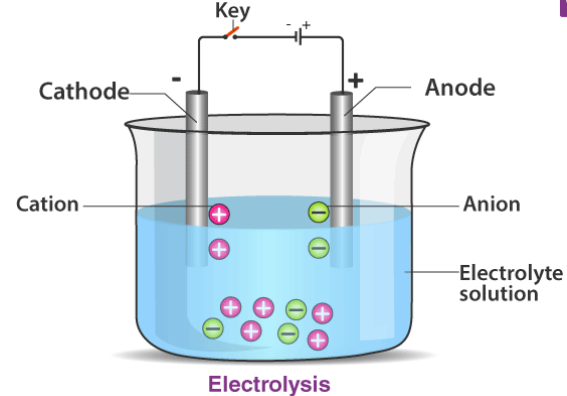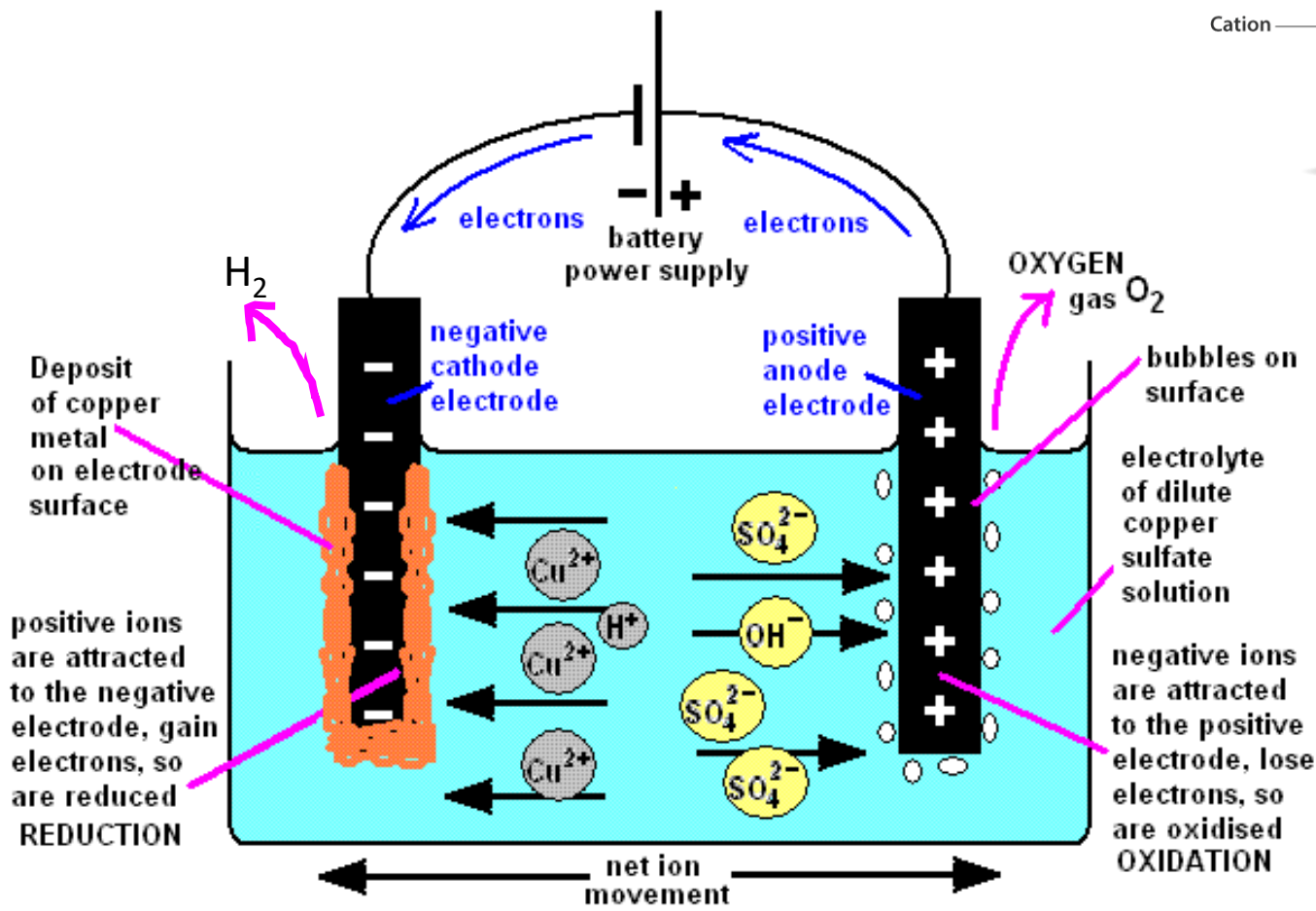

- Cathode - electrode to be coated/plated
- Anode - sacrificial anode or inert anode (platinum or graphite)
- Electrolyte
- DC power supply



# Reaction involved

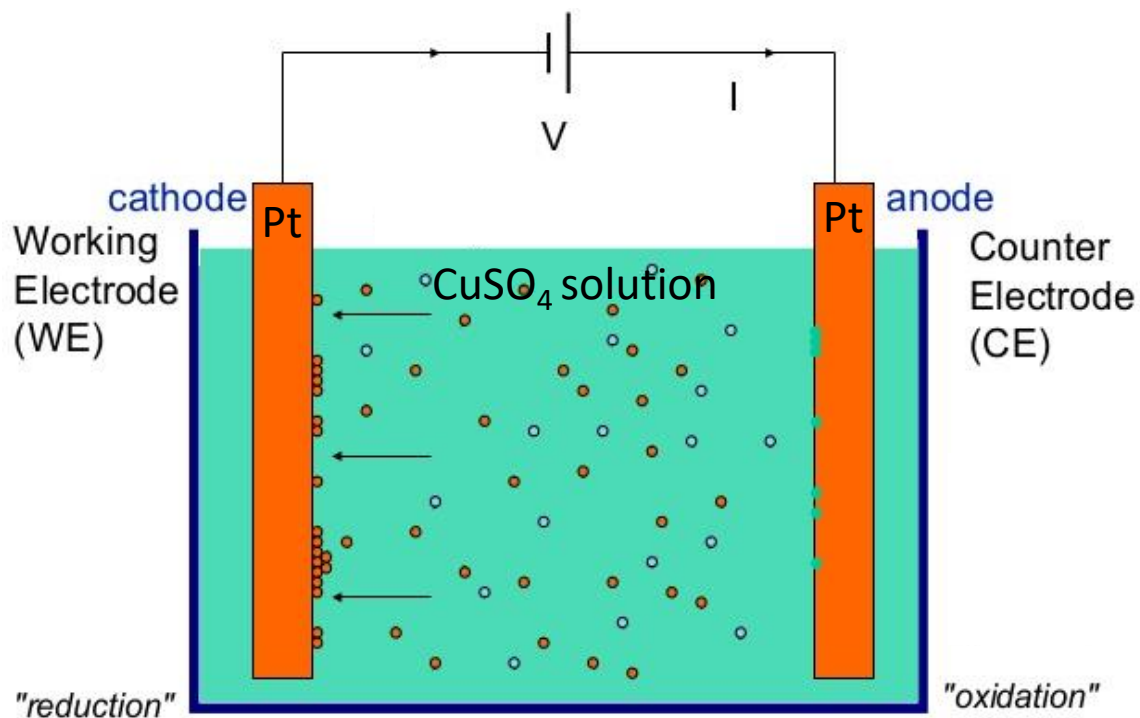


**At Cathode**



**At Anode**

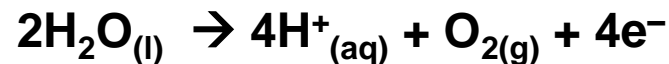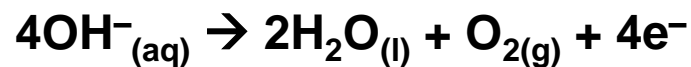

## Circuit Diagram

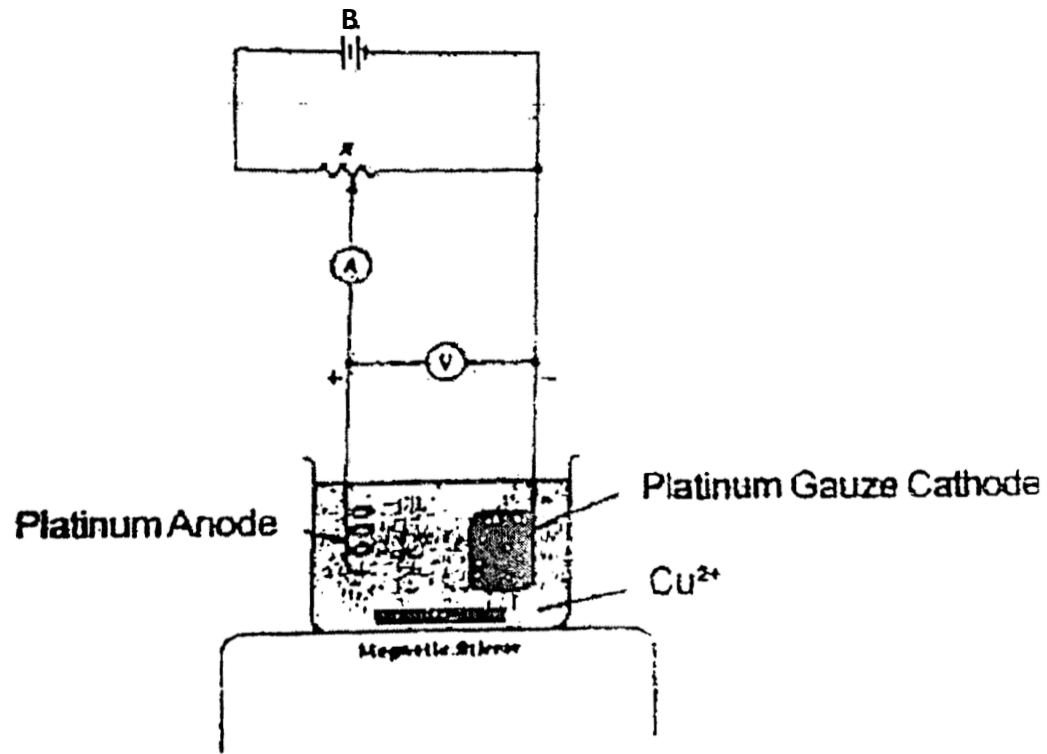

B = Battery.

R = Rheostat.

A = Ammeter,

## PROCEDURE:

- (1) Check the gauge electrode given to you. Consult the instructor if it needs any pretreatment. If necessary give the treatment as it is to be given for spiral platinum anode. Dry it in oven & weight accurately.
- (2) Clean the platinum spiral in 1:1  $\text{HNO}_3$  solution for 3-4 minutes. Wash with tap water and then distilled water.
- (3) Take exactly 150 mL of 2% copper sulphate solution with a burette into 250 mL beaker. The given solution already contains 5 mL conc.  $\text{H}_2\text{SO}_4$  & conc.  $\text{HNO}_3$  which acts as a depolarizer. Add about 50 mL of distilled water.
- (4) Arrange the circuit as shown in the circuit diagram. Be sure that the gauge cathode is connected to the negative terminal & platinum spiral anode is connected to the positive terminal of the power source, as in the circuit diagram. Complete the circuit by plugging in the switch key & raise the beaker containing the electrolyte unit the cathode is completely immersed in the electrolyte. Adjust the rheostat such that the current reading, on the ammeter is about 0.2 ampere. Be sure that the electrodes are not short-circuited & that the magnetic paddle does not touch either of the electrodes.

- (5) Stir the solution vigorously with magnetic stirrer. Continue the electrolysis until blue color of solution has entirely disappeared (This will take about 45 minutes). Reduce the current to 0.1 amp. By adjusting rheostat, add 1.0 g of Urea and continue the electrolysis for another 5 minutes. Test the completion of copper deposition by taking a drop of test solution with glass rod on a filter paper and placing drop of concentrated ammonia solution close to it. Where the boundaries of the two solutions meet, a blue color will be formed if the solution contains copper, if no such blue color is observed; it is an indication that all  $\text{Cu}^{2+}$  ions reduced to Cu and deposited on the cathode.
- (6) To stop electrolysis, turn off magnetic stirrer. Remove the support under beaker and slowly lower the beaker with one hand, while washing the exposed portion of the cathode with distilled water. As soon as the cathode is completely out of solution, cut off the current. Remove the cathode, wash it thoroughly with distilled water and then dip it in a beaker of acetone. Place it on a watch glass and keep it in electric oven for 2 to 3 minutes. Cool the electrode to room temperature in desiccators and then weight accurately

### OBSERVATION:

- (1) Weight of Pt gauze electrode ( $W_1$ ) 12.24 g.  
(2) Weight of Pt gauze + copper deposited ( $W_2$ ) 12.68  
(3) Weight of copper deposited ( $W_3$ ) = ( $W_2 - W_1$ ) = 0.44 g.  
(4) Current employed = \_\_\_\_\_ amp.  
(5) Duration of electrolysis = \_\_\_\_\_ seconds.  
(6) Volume of stock copper solution taken 150 mL.

### CALCULATION:

Copper present in 150 ml of the given stock solution  $W_3$  = 0.44 g.  
Weight of copper per liter in the given stock solution of copper sulphate

$$W_3 \times \frac{1000}{150} = \underline{2.93} \text{ g}$$

### RESULTS:

- (1) 150 mL of the given stock solution contains 0.44 g of Cu.  
(2) Amount of Cu in one liter of the given stock solution = 2.93 g



### NOTES:

1. Acid concentration should not be too high otherwise copper deposition may be incomplete & the deposition may be on adherent.
2. Addition of nitrate ion is necessary for its polarizing action. It acts at cathode as follows  
$$\text{NO}_3^- + 10\text{H}^+ + 8\text{e}^- \rightarrow \text{NH}_4^+ + 3\text{H}_2\text{O}$$
3. Reduction potential of nitrate ion is lower than the discharge potential of  $\text{H}_2$  & hence  $\text{H}_2$  is not liberated.
4. Nitrite ion is formed as follows :  $2\text{H}^+ + \text{NO}_3^- + 2\text{e}^- \rightarrow \text{H}_2\text{O} + \text{NO}_2^-$
5. Urea is added for removal of nitrate, which prevents complete deposition. Nitrite is removed accordingly to  $2\text{NO}_2 + 2\text{H}^+ + \text{CO}(\text{NH}_2)_2 \rightarrow \text{CO}_2 + 2\text{N}_2 + 3\text{H}_2\text{O}$  eqn. Solution should be free from Ag, Bi, Hg, Se, Te, As, Sn, Au, & Pt metals; CNS,  $\text{Cl}^-$  & oxidizing agents, etc